# Einmal-Neutralelektroden / Single-use return plates

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Seite 3

**DEUTSCH**

**Operating Manual**

Seite 11

**ENGLISH**

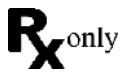
# Inhalt

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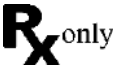
# Contents

1	Symbol descriptions .....	12
2	Introduction .....	13
3	Scope .....	13
3.1	Intended Use .....	13
3.2	Contraindication .....	13
4	Handling .....	13
5	Commissioning the unit.....	14
5.1	Neutral electrode settings.....	14
5.2	Application information .....	14
5.3	Areas for attachment of the neutral electrode .....	15
5.4	Connection to HF unit in monopolar mode.....	16
6	Warnings .....	16
7	Safety Remarks.....	17
8	Protective measures and warnings.....	17
9	Preparation.....	18
10	Packaging .....	18
11	Storage.....	18
12	Warranty/Repair .....	18
13	Servicing and manufacturer address .....	18

# 1 Symbol descriptions

Symbol	Definition
	CE-labelling
	Attention
	Manufacturer
	Lot-description
	Reference code
	Medical device / FDA Prescription device
	Medical device
	Non sterile
	Keep away from sunlight
	Dry storage required
	Temperature
	(Electronic) instruction for use
	No intervention on patients with cardiac pacemakers.

## 2 Introduction

Our products are exclusively intended for professional use by appropriately trained and qualified personnel and may only be acquired by them.

You are now the owner of a high-quality product whose use and correct handling are described in the following. In order to minimize possible risks to patients and users, please observe these instructions carefully. Use, disinfection, cleaning and sterilization may only be performed by suitably trained specialist personnel.

## 3 Scope

The use of HF neutral electrodes takes place in HF surgery.

### 3.1 Intended Use

HF Single use neutral electrode: A non-active conductor attached to a patient and connected to a high frequency electrosurgical generator to establish a circuit for the return of electrical current to the generator following the transmission of the current to perform electrosurgical procedures on the patient. It is attached to the patient where the largest surface can be covered as close to the surgical field as possible. Typically, it has permanently attached cables. It is a disposable item.

### 3.2 Contraindication

The HF neutral electrodes must not be used on damaged and injured skin.

## 4 Handling

The type of treatment must be determined in each individual case by the surgeon in cooperation with the internist and the anaesthetist.

For operational use in various surgical disciplines must be done by appropriately trained and qualified personnel.

## 5 Commissioning the unit

### 5.1 Neutral electrode settings

In the monopolar operating mode, a neutral electrode is required.

After attaching the neutral electrode to the patient, the HF unit recognizes the neutral electrode automatically (cf. illustration on the left and chapter **Fehler! Verweisquelle konnte nicht gefunden werden. Fehler! Verweisquelle konnte nicht gefunden werden.**). The HBS x-touch 200i recognizes neutral electrodes with divided and undivided contact surfaces.

**If the displayed neutral electrode type does not agree with the used neutral electrode type, this is indicative of poor adhesion of the neutral electrode on the patient, a defective neutral electrode or a defective cable!**


**Please note the following information when attaching the neutral electrode**

### 5.2 Application information


- When using single-use products, always pay attention to the use-by date. Do not use electrodes which are over this date!
- Neutral electrodes designed for single use must not be used a second time!
- The electrode must not be attached to surfaces with a high amount of hair growth. If necessary, shave the relevant area before attaching.
- The neutral electrode must be attached in such a way that its entire surface adheres to the patient's body.
- After every change of the patient's position, check that the neutral electrode is adhering correctly and also check the cable routing.
- As soon as the electrode has been attached, it must not be pulled off and attached again. Should the electrode have to be repositioned, use a new neutral electrode.
- The neutral electrode must not come into contact with fluids and must not be coiled round.
- Do not apply any additional electrically conductive gel to the neutral electrode.
- When performing surgery on small children, special electrodes suitable for children or babies must be used .
- Never under any circumstances use damaged products!
- When removing the neutral electrode, never pull on the cable!
- When pulling off single-use electrodes, ensure that the patient's skin is not damaged. Avoid any abrupt movements.

### 5.3 Areas for attachment of the neutral electrode


- Attach the neutral electrode to clean, dry skin.
- The electrode must be attached close to the operation area, but with a minimum direct distance from the operating field of 20 cm.
- Attach over tissue with a good supply of blood without dips or curves in the skin, for instance on the upper arm or thigh. (cf. picture below)
- Do not attach the neutral electrode to projecting surfaces.
- Do not attach above bones, scars, cuts or scratches.
- The electrode must not be attached to surfaces with a high amount of hair growth. If necessary, shave the relevant area before attaching.
- Do not attach above particularly adipose areas such as the abdomen or buttocks.
- The neutral electrode must not be attached above implants.
- Electrodes of monitoring systems must be placed at a distance of at least 20 cm from the operation area and from the position of the attached neutral electrode
- Correct attachment of the neutral electrode in adult patients

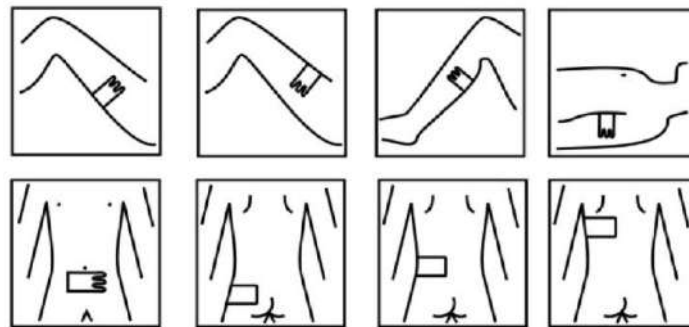- Correct attachment of the neutral electrode in children

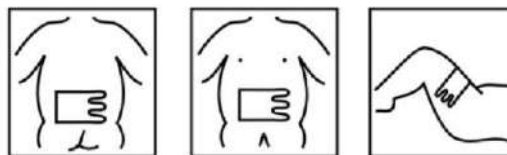- Attachment in patients with pacemakers  
**Note! Consult the treating cardiologist!**

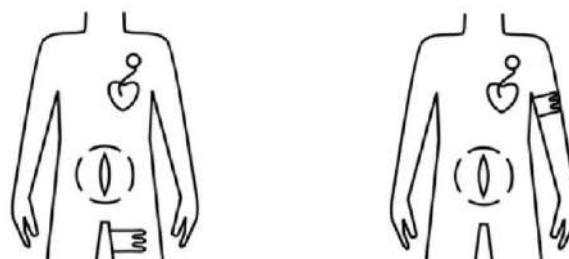**RIGHT**

**WRONG!**

## 5.4 Connection to HF unit in monopolar mode

First fix the neutral electrode ("return plate") to the patient, preferably on the upper arm or on the thigh. Ensure that the skin is free from hair and grease. If a reusable neutral electrode (return plate) is used, do not apply conductive gel. The current paths inside the patient's body should be short and must run diagonally. Current paths must never run in a transverse direction through the body or across the thorax. The patient must be insulated against all electrically conductive parts. Ground the operating table. Place the patient on a dry electrically insulating layer. Prevent skin areas touching each other. Place dry gauze in between. Then connect the neutral electrode, the foot switch and the "active" HF cable to the HF unit.

## 6 Warnings

The products may only be used in combination with HEBU accessories and only by clinically instructed and qualified personnel. Electrodes must not be twisted, as this may cause damage to the material and the electrode might break. Power output of the HF unit: Always select the lowest possible output of the HF unit for any procedure. If the coagulation capacity of the electrode is below normal, do not increase the output of the HF unit without a thorough preliminary check. In particular, the following must be checked: correct attachment of all HF cables and plugs, correct activation of current by the finger - or footswitches, undamaged insulation of the HF cables and instrument, and check that the distal end of the electrode is dean and undamaged. Before operating, check that the signal transmission of the monitors works without interference (e.g. without noise); this check is made by pressing the yellow button (CUT) and the blue button (COAG).


Always unplug the HF-cable by holding on the connector. **Never pull the cable itself**, as this may lead to defects. Cables with defective/brittle insulation or with a broken electrical wire (due to severe bending) might cause a fire. Active HF instruments and HF handles must not be placed on the patient in order to avoid injuries in case of involuntary HF activation or activated "Auto Start" mode.

**Endogenous burns** are those caused by high current density in the patient's tissue. Possible reasons are: The patient is inadvertently positioned in contact with electrically conductive parts. Direct contact between skin areas and HF cables may lead to electrical capacitance which in turn may cause burns.

**Exogenous burns** are those caused by heat of ignited fluids or gases. They may also be caused by explosions. Possible reasons include: Ignition of skin cleaning agents and disinfectants, ignition or narcotic gases etc.

 **Cardiac pacemakers** may be damaged by electrosurgical current. Consult a cardiologist before operating. Never use electrosurgical current on outpatients fitted with pacemakers.
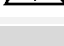
## 7 Safety Remarks

	The use of faulty instruments is in principle forbidden.
	Observe manufacturer's instructions for electrical load capacity!
	Avoid the simultaneous use of HF surgery and electrically grounded electrodes and sensors of diagnostic devices!
	Please observe the additional information enclosed with the products.
	Remove all protective sleeves and films prior to first using or preparation for use.
	The safe combination of different products or of products with implants must be reviewed prior to clinical application by the user.
	Avoid improper throwing or dropping of instruments
	To avoid any damage, HF neutral electrodes with a damaged surface must be removed immediately!
	Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

## 8 Protective measures and warnings

	Use the greatest care when working with the neutral conductor!
	Prevent the penetration of fluids between the patient's body and the neutral electrode.
	A marked drop in power at a normal setting can be an indication of incorrect application of the neutral electrode.

## 9 Preparation

The product is a disposable product and must not be reused.

## 10 Packaging

Compliant packaging of products for sterilization in line with ISO 11607. Packaging used must be suitable for the instruments and protect them from microbiological contamination during storage. The seal must not be under tension. For disposable products, always observe the expiry date. Do not use the electrodes after this date!

## 11 Storage

Storage of the sterile products in a dry, clean and dust-free environment. Products must be protected from direct sunlight in storage.

The storage and transport temperature is 0-40°C.

## 12 Warranty/Repair

Our products are manufactured from high-grade materials and carefully checked prior to dispatch. Defective products should no longer be used.

## 13 Servicing and manufacturer address

Should you require the instructions for use in paper form, please use the contact details below. The instructions for use in paper form will be made available to you within seven calendar days of receipt of the request.

Alternatively, you can print out the electronic instructions for use yourself.

HEBUmedical GmbH  
Badstraße 8  
78532 Tuttlingen / Germany  
Tel. +49 7461 94 71 - 0  
Fax +49 7461 94 71 - 22  
eMail: [service@HEBUmedical.de](mailto:service@HEBUmedical.de)  
Web: [www.HEBUmedical.de](http://www.HEBUmedical.de)

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**Our vast range of products covers over 10 000 various instruments.  
They can be found in Internet or request our catalog.**

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